




# SOLUTIONS

## «DEGENERATIVE SHOULDER»

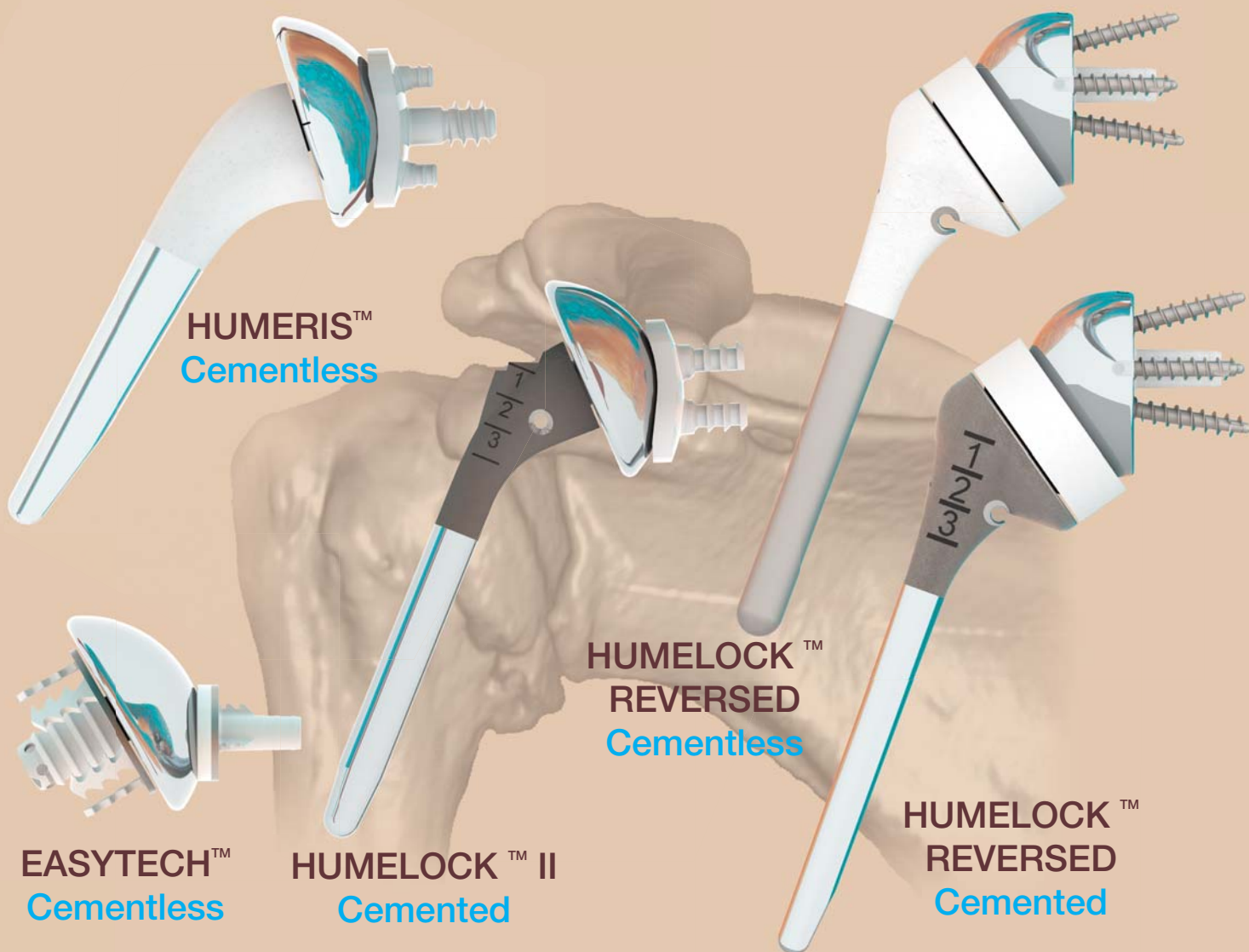

**HUMERIS™**  
Cementless

**HUMELOCK™**  
REVERSED  
Cementless

**EASYTECH™**  
Cementless

**HUMELOCK™ II**  
Cemented

**HUMELOCK™**  
REVERSED  
Cemented

TECHNICAL FLYER

# ALGORITHM




## ANATOMICAL

## REVERSIBLE - REVISION

## REVERSED

Switch of Anatomical Prosthesis to Reversed Prosthesis: remove the humeral head and the double taper, impact an humeral cup 135/145°, transformation of 135° to 145° angle.

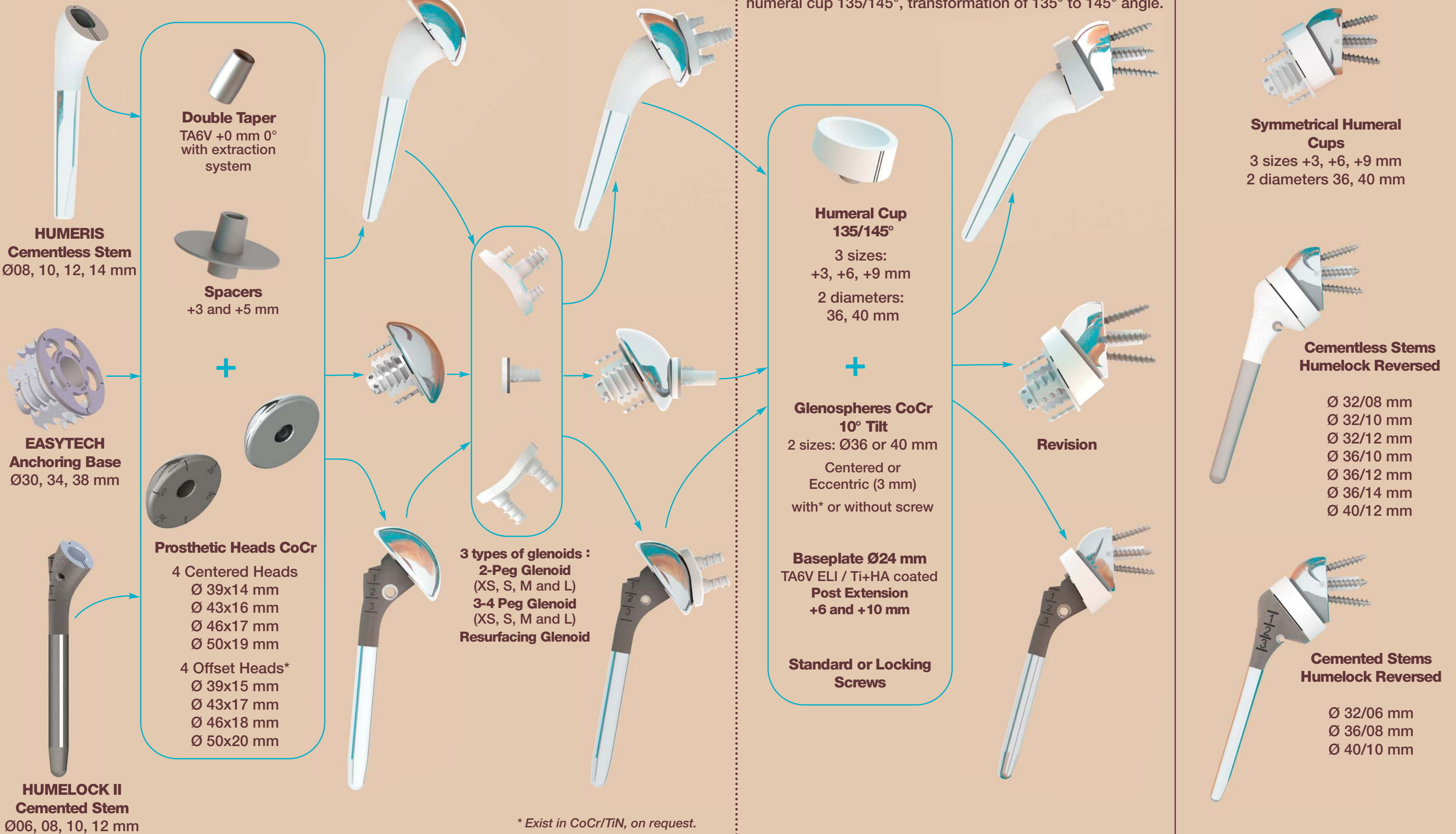

# PROPERTIES

## ANCHORING BASES EASYTECH:

Primary stability from central retentive post and peripheral notches.  
Secondary stability provided by a double layer Ti+HA coating.  
Anti-rotation from alternating notches positioned in different axes / angles.  
Female conical fitting to be reversed.

## Modular System of EASYTECH Prosthesis

<b>Anchoring Base</b>	30	34	38
<b>Head Diam.</b>	39	43	46
	43	46	50



## Humeral Head / Glenoid Association

<b>Glenoid Size</b>	XS	S	M	L
<b>Head Diam.</b>	39	39	43	43
	43	43	46	46
			50	50



## 2-Peg or 3-4 Peg Glenoid:

Made from UHMWPE with radio-opaque marker made from tantalum, 2 retentive pegs with 2 different inter-axis sizes XS-S / M-L. Mismatch included between 5 and 9. Convex, rough base for a better adhesion to cement.

## Cementless Resurfacing Glenoid:

Fits all anatomies, Ø20 mm, net-shape molded UHMWPE inside a TA6V ELI shell, coated with Ti+HA, with 3 anti-rotation fins and extraction system.  
2 post-extensions +6 and +10 mm for longer anchoring.



1663, rue de Majornas - 01440 Viriat - France  
Tel. : (33) 04 74 55 35 55 - Fax : (33) 04 74 52 44 01  
E-mail: info@fxsolutions.fr - www.fxsolutions.fr